


The vast majority of Implantologists & Periodontists recommend PTFE sutures (Profimed class) for the best results:

- ▼ Easy handling
- ▼ Soft & pliable
- ▼ Minimal tissue trauma & drag
- ▼ Resistance to microbes harboring
- ▼ Minimal to zero package memory
- ▼ Biologically inert
- ▼ Excellent knot security
- ▼ No discomfort for the patient (soft edges)
- ▼ Ideal aesthetic results
- ▼ No degradation
- ▼ No weak spots on the strand when crimped
- ▼ No additives, pure PTFE

 **medipac**
dental surgery

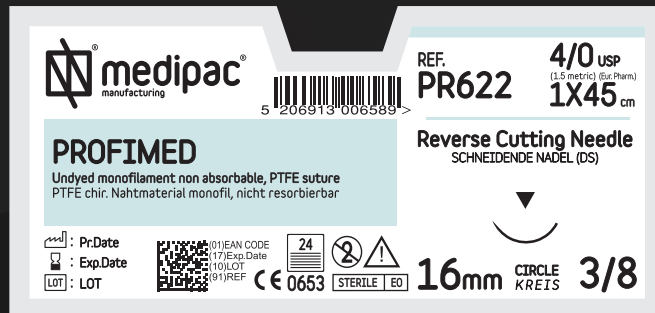

medipac
manufacturing

REF. **PR622** **4/0 USP**
(1.5 metric) (6x Pharm)
1X45 cm

PROFIMED
Undyed monofilament non absorbable, PTFE suture
PTFE chir. Nahtmaterial monofil, nicht resorbierbar

Reverse Cutting Needle
SCHNEIDENDE NADEL (DS)

Pr.Date
Exp.Date
LOT

01)EAN CODE
07)Exp.Date
10)LOT
09)REF

24
STERILE EO

16mm **CIRCLE** **3/8**
KRETS

 **medipac**
dental surgery

Medipac S.A.
Industrial Area P.O. Box 1, 611 00, Kilkis, Greece
T: +30 23410 71991, F: +30 23410 71979
E: info@medipac.gr, W: www.medipac.gr

Are you looking for Results?

PROFIMED
is the Answer

Look what you have been missing out!

PROFIMED

We have meticulously prepared the following Table for you

Profimed is a sterile non-absorbable monofilament surgical suture made from polytetrafluoroethylene (PTFE) polymer. The suture is available undyed in a wide variety of needle-thread combinations.

Suture category	Package Memory	Tissue drag & trauma	Resistance to microbes harbouring - minimal bacterial plate on the strand surface	Possible weak spots on the strand through crushing	Easy handling - Pliability - Softness	Knot security	Degradation - Gradual loss of tensile strength	Probability of creating scars on the skin	Discomfort in the oral cavity caused by the edges
Profimed (PTFE)	No	Minimal	Yes	No	Yes	Yes*	No	Less	No
Polyamid	Yes	Minimal	Yes	Yes	No	Less	Extremely slow	Less	Yes
Propylen	Yes	Minimal	Yes	Yes	No	Less	No	Less	Yes
Silk	No	More	Minimal	No	Yes	More	Yes (Very slow)	More	No
Polyester	No	More	Low	No	Yes	More	No	More	No
PGA	No	More	Low	No	Yes	More	Yes	More	No
Neosorb (PGLA)	No	More	Low	No	Yes	More	Yes	More	No
Neosorb Rapid (RPGLA)	No	More	Low	No	Yes	More	Yes	More	No
Monofast (PGCL)	Yes	Minimal	Yes	Yes	Medium	Less	Yes	Less	Medium
Monosorb (PDO)	Yes	Minimal	Yes	Yes	No	Less	Yes	Less	Yes

*(special technique required)

All Medipac sutures are available in wide variety of thread length, diameter and needle combinations. Please visit our online catalogue at www.medipac.gr